

Reprocessed by



Sustainability Solutions

Instructions for Use **Reprocessed Posey® Chair Alarm Sensor Pad**

Reprocessed Device for Single Patient Use

Caution: Federal (U.S.A.) law restricts this device to sale by or on the order of a physician.

- NON-STERILE
- Exposed to Ethylene Oxide (EO) gas

Explanation of Symbols



Date of Reprocessing



Do Not Reuse



See Instructions For Use

Stryker Sustainability Solutions, Inc. © 2019
1810 W Drake Dr.
Tempe AZ, 85283
sustainability.stryker.com
888.888.3433

Reprocessed Posey® Chair Alarm Sensor Pad Description

The Reprocessed Posey® Chair Alarm Sensor Pads are sensor pads designed for use with Posey® Fall Alarms to detect when a patient has taken weight off the sensor to help notify nurses and caregivers that the patient has gotten up from a sitting or resting position.

Indications for Use

The Stryker Sustainability Solutions Reprocessed Posey® Chair Alarm Sensor Pads alert staff that a patient has taken weight off the sensor and may be at immediate risk of falling. However, fall alarm systems DO NOT prevent falls or injury from falls and are not a substitute for patient care, caregiver rounding and a comprehensive fall management protocol in your facility.

Warnings

- ALWAYS check sensor pad when connecting it to a Posey® alarm. You can check a sensor pad by attaching it to the sensor input on the alarm, activating the alarm and placing pressure on the sensor pad. When the pressure is released, the alarm should sound. Repeat this pressure/release test in several different areas along the entire length of the sensor pad to ensure entire sensor pad functions properly. If the alarm and/or sensor pad do not function properly, remove the alarm and sensor pad from service and replace them with a properly functioning alarm and/or sensor pad. DO NOT use the alarm or sensor pad if it does not activate each time weight is removed from the sensor pad.
 - Note: Seating or positioning aids or heavy cushions may keep alarm from activating. If the test fails, remove or relocate aids, or try a different position for sensor such as on top of cushion.
 - If the sensor fails any test, stop use at once and replace with a new sensor.
 - Sensor pads past warranty expiration date may fail, causing false alarm or no alarm.
- NEVER jerk or pull on cord to remove plug. Doing so will damage the cord or plug and may cause sensor pad to fail. ALWAYS use plastic tab to release plug.
- NEVER roll, fold, or crease sensor pad. This will damage sensor pad and may cause false alarm or no alarm.
- To reduce the risk of serious injury, ALWAYS follow the Directions for Testing Sensor Pad (see below) after putting the sensor pad in place and before leaving patient unattended. DO NOT use any alarm or sensor pad that does not alarm each time it is tested.
- USE OF RESTRAINTS: The position of the sensor pad is vital when using a restraint. Make sure the restraint is applied correctly per instructions for that device. Straps must NOT cross over sensor pad.
 - If straps cross over sensor pad and patient moves, pressure from straps may prevent alarm from activating. If patient falls out of chair and is suspended in the restraint, serious injury may occur from chest compression or suffocation.
 - Use extreme caution with chair cushions. If a cushion dislodges, straps may loosen and allow patient to slide off seat and become suspended.
- The sensor pad will not work if plugged into phone outlet.

Precautions

- This device may not be suitable for all high fall-risk patients.
- To prevent damage, NEVER roll, fold, or crease sensor.
- NEVER connect Posey® sensors to other manufacturers' alarm units.
- Sensor pads beyond expiration date are out of warranty and may not function properly, causing false alarm or no alarm. When putting sensor pad into use, record the first use date and calculate expiration date from initial use.
- To avoid inconvenience to staff and patients, and to protect sensor pads from damage, you should follow these steps:
 - Only use Reprocessed Posey® sensor pads with the Posey® alarm.
 - When routing sensor pad cord to alarm, check that there is slack in the cord.
 - NEVER jerk or pull on the cord to remove plug. Doing so will damage cord wires or plug.
 - ALWAYS use the plastic tab to release plug.
 - Make sure sensor pad intake ("neck" of chair sensor pad) is clear and not blocked. Air must flow freely in and out of sensor pad for alarm to function. Make sure liquid does not enter at "neck" of sensor pad, as this will damage sensor pad. If needed, use an incontinence pad to protect sensor pad from urine or other liquids.
 - NEVER roll, fold or crease sensor pad. This will damage sensor pad and may cause false alarm or no alarm.
 - Make sure sensor pad lays FLAT on chair surface, directly under patient's weight, and that sensor pad cord is not folded back under the sensor pad.

Directions for Use

1. Read all sensor pad and alarm instructions prior to use.
2. Check that sensor pad, cord, and plug are clean and undamaged. If you have a damaged sensor pad and/or plug, the sensor pad unit should be discarded and replaced.
3. When putting sensor pad into use, record the first use date and calculate expiration date from initial use.

4. Place sensor pad FLAT across width of seat, labeled side UP (refer to Figure 1). NEVER roll, fold, or crease sensor pad. Make sure sensor pad cord is either to the back or side of seat and is not folded back under the sensor pad.
 - Note: If a seat cushion is used, place sensor pad ON TOP of cushion.
3. Adjust so sensor pad is directly under patient's buttocks when seated.
 - Sensor pad should be towards FRONT of chair seat if patient normally sits toward front.
 - Sensor pad should be towards BACK of chair seat if using a posture support or if patient is at risk of forward sliding.
4. If needed, use an incontinence pad to prevent damage from urine or other liquids. If liquid enters at "neck" of sensor pad, alarm may not function.
5. Route the sensor pad cord to the Posey® alarm. Check that the sensor pad cord has slack, is not stressed, is clear of moving parts of chair, and does not pose a tripping hazard.
6. Insert plug into sensor input on the alarm. Check that that the plug on the sensor pad cord is not damaged (plug broken, or wires disconnected) and is securely connected to the alarm.
7. Test sensor pad and alarm (refer to Directions for Testing Sensor Pad).

Directions for Testing Sensor Pad

1. Make sure alarm is ON and in monitoring mode.
2. Check that the plug on the sensor pad cord is not damaged (plug broken, or wires disconnected) and is securely connected to the alarm.
3. Disconnecting the sensor pad from the alarm when the power is on will cause the alarm to activate. This is called a "failsafe" mode. Disconnect the sensor pad to make sure the failsafe mode works. DO NOT use the alarm if the alarm does not sound when the sensor pad is disconnected.
4. When connecting the alarm to the nurse call system, check that the nurse call cable is securely connected to the alarm and the nurse call panel. ALWAYS test alarm and nurse call function if nurse call cable is plugged into the alarm and wall jack. Activate the alarm (remove pressure from sensor pad) and make sure the nurse call light for the proper room activates in the appropriate nurse's station location. Remove the cable from the wall jack and make sure the visual or audible alert at the nurse's station immediately activates.
5. Inspect sensor pad cord and nurse call cable (if in use) to ensure they are out of the footpath and DO NOT pose a tripping hazard.
6. If at any time the sensor pad does not function properly during testing prior to each use, stop use immediately and replace with a new sensor pad.

Storage and Handling

- Store sensor pads flat or hang in a dry, secure environment. DO NOT roll, bend or fold sensor pads, as it may damage internal electronic parts and cause a malfunction.
- This device is designed for use in normal indoor environments.
- This device may be stored in ambient warehouse temperatures at normal humidity levels (10 to 50%). Avoid excess moisture or high humidity that may damage product materials (greater than 90%).

Compatibility

- The Reprocessed Posey® Chair Alarm Sensor Pads are designed for CHAIR USE ONLY with Posey® Fall Alarms only. Never connect the Reprocessed Posey® Chair Alarm Sensor Pad to other manufacturers' sensors and alarms.

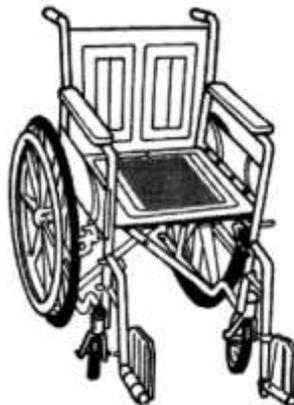


Figure 1 – Posey® Chair Alarm Sensor Pad Positioning

Warranty

8308 – This product is warranted for a period of six (6) months from date of first use.

8309 – This product is warranted for a period of thirty (30) days from date of first use.

8309EL - This product is warranted for a period of thirty (30) days from date of first use.

Reprocessed Products

Unless agreed otherwise in writing, Stryker warrants all reprocessed products, subject to the exceptions provided herein, to be free from defects in reprocessing and to substantially conform to the product specifications contained in the documentation provided by Stryker with the products for one use in accordance with the instructions for use of such product.

STRYKER SHALL NOT BE LAIBLE FOR ANY DAMAGES TO THE EXTENT CAUSED BY ANY DEFECT IN MATERIAL, WORKMANSHIP OR DESIGN BY THE ORIGINAL MANUFACTURER OF THE PRODUCT OR ANY ACT OR OMISSION OF THE ORIGINAL MANUFACTURER OF THE PRODUCT.

Products for which Stryker is the Original Manufacturer

Stryker warrants all products for which it is the original manufacturer, subject to the exceptions provided herein, to be free from defects in design, materials and workmanship and to substantially conform to the product specifications contained in the documentation provided by Stryker with the products for a period of one year from the date of purchase.

General Warranty Terms Applicable to All Products

TO THE FULLEST EXTENT PERMITTED BY LAW, THE EXPRESS WARRANTY SET FORTH HEREIN IS THE ONLY WARRANTY APPLICABLE TO THE PRODUCTS AND IS EXPRESSLY IN LIEU OF ANY OTHER WARRANTY BY STRYKER, EXPRESSED OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, ANY IMPLIED WARRANTY OR MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. IN NO EVENT WILL STRYKER'S LIABILITY ARISING IN CONNECTION WITH THE SALE OF THE PRODUCT (WHETHER UNDER THE THEORIES OF BREACH OF CONTRACT, TORT, MISREPRESENTATION, FRAUD, WARRANTY, NEGLIGENCE, STRICT LIABILITY OR ANY OTHER THEORY OF LAW) EXCEED THE PURCHASE PRICE, CURRENT MARKET VALUE OR RESIDUAL VALUE OF THE PRODUCTS, WHICHEVER IS LESS. STRYKER SHALL NOT BE LIABLE FOR INDIRECT, SPECIAL, INCIDENTAL, PUNITIVE OR CONSEQUENTIAL DAMAGES RESULTING FROM ANY BREACH OF WARRANTY OR UNDER ANY OTHER LEGAL THEORY.

This warranty shall apply only to the original end-user purchaser of products directly from Stryker or a Stryker authorized distributor. This warranty may not be transferred or assigned without the express written consent of Stryker.

This warranty does not apply to: (1) products that have been misused, neglected, modified, altered, adjusted, tampered with, improperly installed or refurbished; (2) products that have been repaired by any person other than Stryker personnel without the prior written consent of Stryker; (3) products that have been subjected to unusual stress or have not been maintained in accordance with the instructions in the user manual or as demonstrated by a Stryker representative; (4) products on which any original serial numbers or other identification marks have been removed or destroyed; (5) products that have been repaired with any unauthorized or non-Stryker components; or (6) failure to follow the instructions for use.

If a valid warranty claim is received within thirty (30) days of the expiration of the applicable warranty period, Stryker will, in its sole discretion: (1) replace the product at no charge with a product that is at least functionally equivalent to the original product or (2) refund the purchase price of the product. If a refund is provided by Stryker, the product for which the refund is provided must be returned to Stryker and will become Stryker's property. In any event, Stryker's liability for breach of warranty shall be limited to the replacement value of the defective or non-conforming part or component.

If Stryker determines in its reasonable discretion that the claimed defect or non-conformance in the product is excluded from warranty coverage as described hereunder, it will notify the customer of such determination and will provide an estimate of the cost of repair of the product. In such an event, any repair would be performed at Stryker's standard rates.

Products and product components repaired or replaced under this warranty continue to be warranted as described herein during the initial applicable warranty period or, if the initial warranty period has expired by the time the product is repaired or replaced, for thirty (30) days after delivery of the repaired or replaced product. When a product or component is replaced, the item provided in replacement will be the customer's property and the replaced item will be Stryker's property. If a refund is provided by Stryker, the product for which the refund is provided must be returned to Stryker and will become Stryker's

property.

The OEM information listed on the label is provided as device ID prior to reprocessing and may contain the trademarks of unrelated third parties that do not sponsor this device.

This product and its packaging have been exposed to ethylene oxide gas (EtO). Even though the product then is processed in compliance with all applicable laws and regulations relating to EtO exposure, Proposition 65, a State of California voter initiative, requires the following notice:

Warning: This product and its packaging have been exposed to ethylene oxide. The packaging may expose you to ethylene oxide, a chemical known to the State of California to cause cancer or birth defects or other reproductive harm.

Posey® is a registered trademark of Posey Company.

HAS EL10103 Rev. B 09/2019