

Gynnie®

OB/GYN Stretcher




Easy. Multi-functional. Convenient.

Gynnie combines the durability and features of Stryker's standard stretcher with supreme OB/GYN functionality. This multi-functional stretcher is ideal for pelvic exams and can also be utilized for general use.



Pelvic exam conversion

3 easy steps to convert stretcher into a pelvic exam with self-locking stirrups.

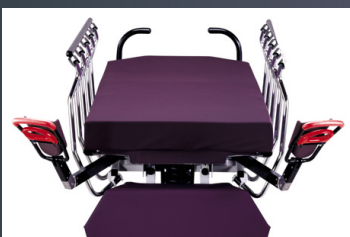

4" Ultra Comfort support surface

Heat sealed seams help prevent fluid collection.



Glideaway siderails

Swings up and out for patient comfort and allows for zero transfer gap when lowered.

Retractable fifth wheel helps improve mobility.

Permanently attached parts to maximize caregivers' time.

Optional calf supports

Dual pedestal hydraulics to help optimize stability.

Four wheel, steel ring brakes with dual side and head activators to help secure platform and accessibility.

Gynnie

OB/GYN Stretcher



Standard features

- 4" (10 cm) Ultra Comfort support surface with sealed seams
- 8" (20 cm) Omni Surface casters with wheel covers
- 26" (66 cm) wide patient surface
- Dual cylinder pneumatic backrest
- Dual pedestal hydraulics with constant rate descent
- Dual side-mounted foot control hydraulics
- Retractable fifth-wheel steering
- Four-wheel, steel-ring brakes with dual side and head activators
- Glideaway siderails
- Integrated storage tray with oxygen bottle holder
- Permanently-attached folding foot section
- Permanently-attached foot supports
- Pop-up push handles
- Roller bumpers
- Two IV receptacles
- Unique one-handed foot section release lever

Stryker reserves the right to change specifications without notice.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Gynnie, Ultra Comfort, Stryker. All other trademarks are trademarks of their respective owners or holder.

3800 E. Centre Avenue
Portage, MI 49002 USA
t: 269 329 2100
toll free: 800 327 0770

patienthandling.stryker.com

Specifications

Model number	1061
Overall length	81" (205.7 cm)
Overall width	
Siderails up	31.5" (80 cm)
Siderails down	28" (71 cm)
Weight capacity	500 lb. (225 kg)
Height range (to litter top)	
High	35.75" (90.8 cm)
Low	22" (55.5 cm)
Litter positioning	
Backrest	0°– 90°
Knee gatch	0°– 90°
Trend./Reverse Trend.	±18°
Maximum angle	87.5°
Patient surface	26" x 75.5" (66 x 191.5 cm)
Siderails	13" x 39" (34 x 99 cm)
Caster diameter	8" (20 cm)
Break away point from wall	68" (172.7 cm)

Optional features

- Hoop for drainage basin
- Patient restraints
- Serving/instrument tray
- Siderail pads
- IV poles
 - Permanent folding two- or three-stage at foot end
 - Removable
- Upright oxygen bottle holder
- Extended warranty available intermediate supports
- Upright oxygen bottle holder

Warranty

One year parts, labor and travel or two years parts