

Symmetry Plus

Treatment Recliner



Recliners are an essential part of every patient room, giving patients freedom from their bed. Our Treatment Recliner features an infinite number of recline positions and full flat-bed position.



Pull-out foot tray gives patients comfort while in transit.

Stowable 360 degree swivel sidetray option.

Wide arm caps are extended to facilitate easy entry and exit.

Wide or narrow widths are available.

Rigid push bar for easy manoeuvrability.


226.8kg (500 lbs) weight capacity warranted.

Backrest locks into infinite number of recline positions.


Full flat-bed position maximises patient comfort.

BackSmart Trendelenburg positioning.

Wide seat comfortably accommodates larger patients.

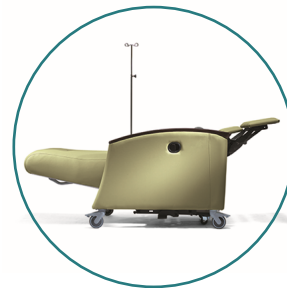

Trendelenburg position

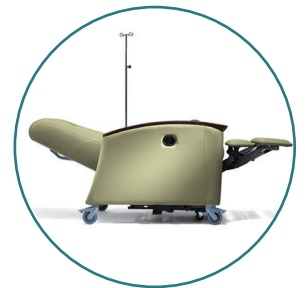

Full flat-bed position

Step-lock brakes



4-sided brake access



Optional side tray



Symmetry Plus Treatment Recliner

Key features

- Independent ottoman allows feet to be raised without moving the back
- Mid-ottoman cushion provides extra support when feet are elevated
- Chair construction includes metal seat frame and sinuous wire springs
- Highly-resilient polyurethane foam cushions
- Removable back cushion enables thorough cleaning and service
- Pull-out foot tray gives patients added comfort while in transit
- Rigid push bar for easy manoeuvrability
- Fixed sides provide stability for electromechanical mounts
- Egress assist polyurethane arm caps
- Meets flammability standards of BIFMA F1-1978 - revised 1980 and CAL 117
- Meets ANSI/BIFMA X5.4 2005 per tests conducted in a certified test facility
- Oversized tracking, lockable casters provide agile manoeuvrability (front rigid, rear swivel)
- **Limited lifetime warranty**

Optional features

- Fold-down side tray provides prep area for caregivers and serves patients' ancillary needs (permanently mounted to right side, left side or both)
- Rear-mounted IV pole and bracket save valuable floor space and allow for easier transport
- Egress assist wood arm caps
- Hinged push bar for manoeuvrability
- Four-wheel swivel with locking casters
- Choose from a variety of stains and fabrics
- CAL 133
- Foley bag hooks on patient left, patient right, or both sides.

* BackSmart Trendelenburg pedals add 5 cm to overall width of chair for clearance.

Specifications

Model	Standard width 3500-000-710	Narrow width 3500-000-705
Overall width	83.8 cm (33")	76.2 cm (30")
Overall height	119.4 cm (47")	119.4 cm (47")
Depth		
Overall	96.5 cm (38")	96.5 cm (38")
Seat	52.1 cm (20.5")	52.1 cm (20.5")
Recline/legrest up	188.0 cm (74")	188.0 cm (74")
Weight	74.8 kg (165 lb.)	77.1 kg (170 lb.)
Seat width	66.0 cm (26")	58.4 cm (23")
Back height from seat	78.7 cm (31")	78.7 cm (31")
Seat height from floor	52.1 cm (20.5")	52.1 cm (20.5")
Arm height from floor	71.1 cm (28")	71.1 cm (28")

A healthcare professional must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that healthcare professionals be trained in the use of any particular product before using it. The information presented is intended to demonstrate the breadth of Stryker product offerings. A healthcare professional must always refer to the package insert, product label and/or instructions for use before using any Stryker product. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area. Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Stryker. All other trademarks are trademarks of their respective owners or holders. For more information please contact your local Stryker sales representative.