

# T2 Alpha

Nailing redefined.





**Stability** when you need it

# T2 Alpha

## Nailing redefined.

T2 Alpha introduces a new screw option called the Advanced Locking Screw designed for added axial and angular stability. Advanced Locking Screws are designed with oversized threads that engage with internal threads of the T2 Alpha nails while maintaining bicortical purchase.<sup>1</sup> With the ability for use in any round 5mm hole of T2 Alpha nails, Advanced Locking Screws allow surgeons to choose when additional stability is needed.

- ▶ Advanced Locking Screws thread and lock into any round 5mm hole
- ▶ Self-centering screw tip designed for simplified insertion
  - 18° cone of angulation for insertion
- ▶ Designed to limit relative movement between nail and screw<sup>2</sup>

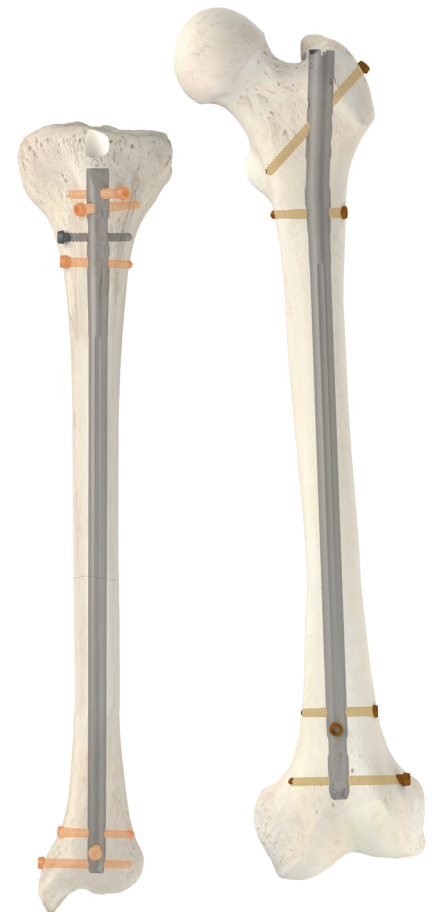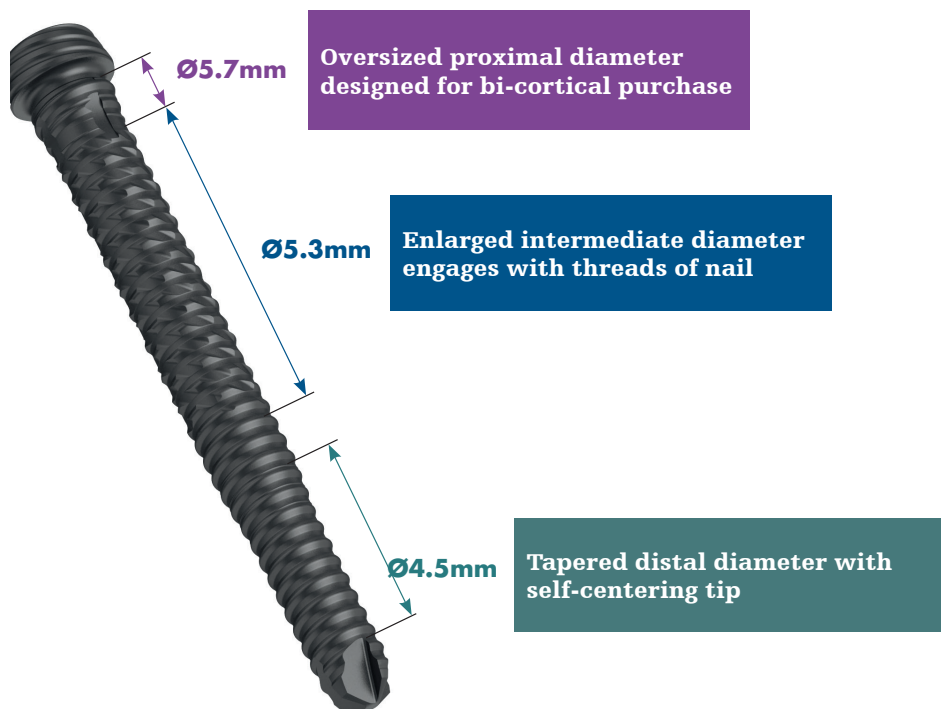

Colors indicate where Advanced Locking Screws can be utilized. Advanced Locking Screws are not accepted in the oblong or recon holes of the nail.

### References:

1. T2-WP-6, 07-2018
2. T2 Alpha Femur Antegrade GT/PF Nailing System Operative Technique, T2-ST-17\_Rev-4, 09-2018

## Trauma & Extremities

This document is intended solely for the use of healthcare professionals. A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. We do not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate Stryker's products. A surgeon must always refer to the package insert, product label and/or instructions for use, including the instructions for cleaning and sterilization (if applicable), before using any of Stryker's products. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your representative if you have questions about the availability of Stryker's products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Stryker, T2 Alpha. All other trademarks are trademarks of their respective owners or holders.

Content ID: T2-SS-7, 01-2019  
Copyright © 2019 Stryker

**Stryker Australia**  
8 Herbert Street,  
St Leonards NSW 2065  
T: 61 2 9467 1000 F: 61 2 9467 1010

**Stryker New Zealand**  
PO Box 17-136  
Greenlane Auckland 1546  
T: 64 9 573 1890 F: 64 9 573 1891

[stryker.com](http://stryker.com)