


# DartFire Edge<sup>®</sup>

## Cannulated Screw System

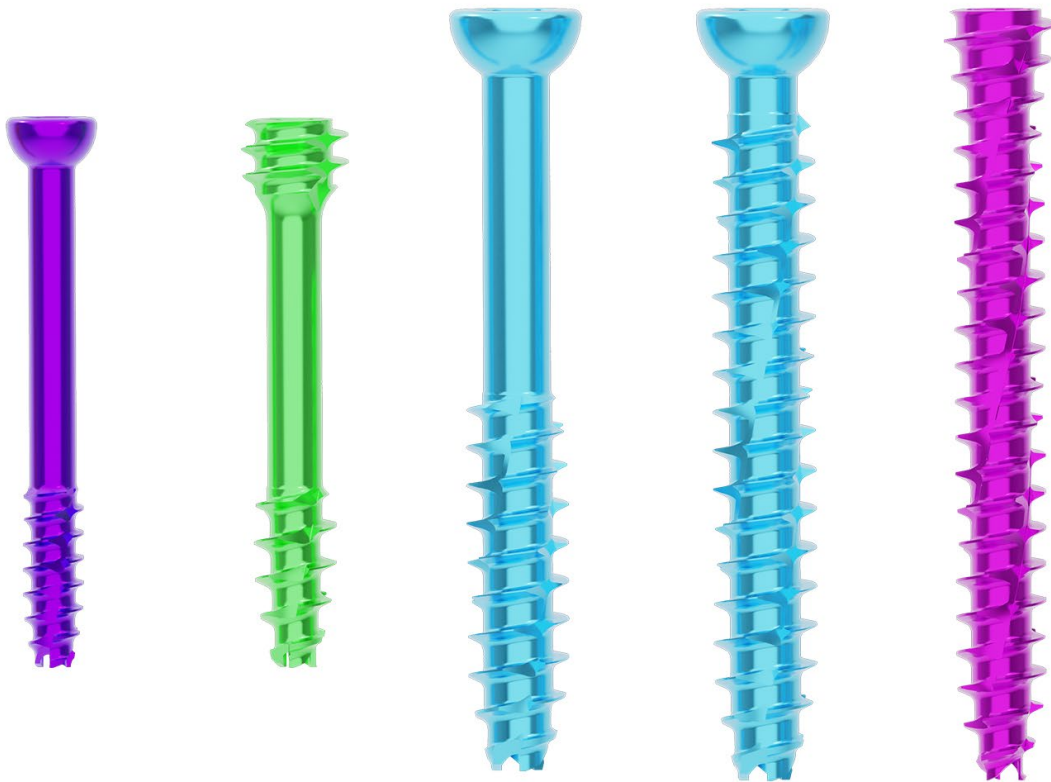

# DartFire Edge®

## Cannulated Screw System

The **DartFire Edge Cannulated Screw System** is a comprehensive forefoot and midfoot screw offering delivered in a ready-for-surgery configuration. Don't settle for what's available. Stop dealing with SPD surprises and delays. Choose a universal, sterile packaged system that fits your practice—a whole new routine driven by efficiency and personalization.

**DartFire Edge** was designed to include all the features of multiple systems available to you in one universal offering. It's goal, to simplify inventory management by only needing to stock a single screw system with everything you need. Intended to save time and space, you no longer need multiple trays of screws in the OR. The simplified and sterile packaging design allows your staff to no longer worry about processing each tray through the sterilization department. Grab the screw needed along with its associated instrument pack and you are ready for the case.

### DartFire Edge provides

- Low profile, self-drilling and self-tapping, titanium alloy (Ti 6Al-4V) compression screws designed to provide exceptional bite and reliable sharpness with cutting flutes running the length of the thread
- Reverse cutting flutes designed to aid in ease of removal
- Headed and headless options ranging from 2.0mm to 4.0mm
- Short thread, long thread, full thread options
- Sterile instrument pack with K-wires, cannulated drill, cannulated depthsink, driver, handle (AO connection)
- 2.0mm and 2.7mm snap-off screws
- Color-coded implant and instrument pack labeling for easy identification



## Trayless efficiency

DartFire Edge takes away the headache of a cluttered OR table from instruments and disposables across multiple trays

- All screws and instrument packs come individually packaged and pre-sterilized
- Easy-open tabs
- See-through windows
- High visibility, color-coded labeling
- 72% decreased footprint compared with existing packaging<sup>1</sup>
- 82% reduction in OR footprint compared to reusable trays<sup>2</sup>
- Outer packaging layers are 100% recyclable



## Simply sharp

A screw system reimagined, not reprocessed

- Single use screws and instruments provide reliable sharpness for every case
- Cutting flutes run the entire length of the thread creating an exceptional bite



## Designed for your practice

Only stock what you need

- As part of the **Pulse Intelligent Delivery Platform**, procedure mapping technology allows offered diameters and lengths of screws to be personalized to your usage
- Always have what you need and nothing that you don't
- A universal screw offering with headed and headless options ranging from short, long, and full thread giving your practice one reliable system

# Most comprehensive cannulated screw system in **Stryker's Foot & Ankle portfolio**

## Headed options

Type	Diameter	Length
Short thread	2.0mm	10-40mm
Short thread	2.5mm	10-40mm
Short thread	3.0mm	12-50mm
Short thread	3.5mm	12-50mm
Short thread	4.0mm	10-40mm

Long thread	4.0mm	16-50mm
-------------	-------	---------

Full thread	3.5mm	12-50mm
Full thread	4.0mm	12-50mm

## Headless options

Type	Diameter	Length
Short thread	2.0mm	10-40mm
Short thread	2.5mm	10-40mm
Short thread	3.0mm	12-50mm
Short thread	3.5mm	12-50mm
Short thread	4.0mm	10-40mm

Long thread	4.0mm	16-50mm
-------------	-------	---------

Full thread	2.5mm	10-40mm
Full thread	3.0mm	12-50mm
Full thread	3.5mm	12-50mm
Full thread	4.0mm	12-50mm

## Snap-off options

Diameter	Length
2.0mm	11mm, 12mm, 13mm
2.7mm	13mm, 15mm, 17mm

\*Each length increases in 2mm increments

## Foot & Ankle

Visit [footankle.stryker.com](https://footankle.stryker.com) or email us at [footandankle@stryker.com](mailto:footandankle@stryker.com)

1. Data on file
2. Space and cost savings may vary. Data on file on ASC Case study. For further information, reference AP-017256.

The information contained in this document is intended for healthcare professionals only. A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker's product offerings. A surgeon must always refer to the package insert, product label and/or instructions for use before using any of Stryker's products. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your sales representative if you have questions about the availability of products in your area. Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: DartFire Edge, Stryker. All other trademarks are trademarks of their respective owners or holders.