


## Superior strength\* and profile for uncompromised performance

### Strength\*

Type II Ti alloy provides increased material fatigue strength superior to color anodized Ti alloy and stainless steel

### Profile


Low-profile plates designed to reduce soft tissue irritation

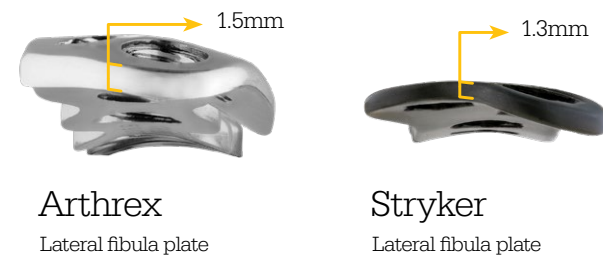

14% Decrease in plate thickness vs Arthrex lateral fibula plate\*\*

### Speed

Implants and instruments designed to reduce steps and improve efficiency

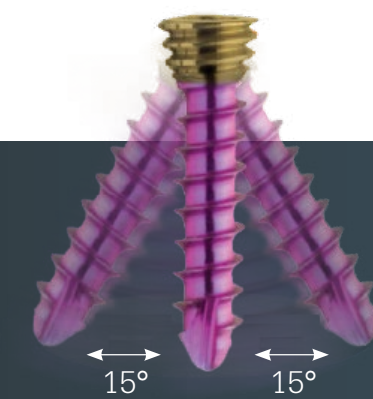


## 3Di technology

### Unsurpassed polyaxial locking

- Titanium Nitride coated heads coupled with proprietary plate hole design strikes a perfect balance for an unsurpassed locking mechanism
- Off axis locking up to 15°, in any direction, with the ability to relock screws up to 3 times
- Available in 2.7mm and 3.5mm diameters

\* Data on file.  
\*\* Based on publicly available product literature

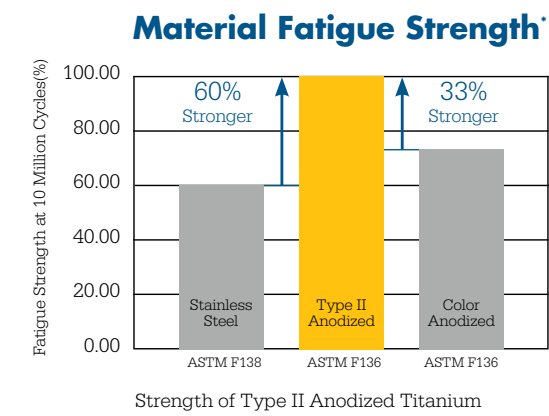

# Ortholoc® 2

with 3Di technology

Superior strength\* and profile for uncompromised performance



### Foot & Ankle

This document is intended solely for the use of healthcare professionals. A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. We do not dispense medical advice and recommend that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate Stryker's products. A surgeon must always refer to the package insert, product label and/or instructions for use, including the instructions for cleaning and sterilization (if applicable), before using any of Stryker's products. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your representative if you have questions about the availability of Stryker's products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Ortholoc, Stryker. All other trademarks are trademarks of their respective owners or holders.

The absence of a product, feature, or service name, or logo from this list does not constitute a waiver of Stryker's trademark or other intellectual property rights concerning that name or logo.

Content ID: 016405A  
Copyright © 2021 Stryker



**Manufactured by:**  
Wright Medical Technology, Inc.  
1023 Cherry Road  
Memphis, TN 38117

stryker.com

# Ortholoc 2

with 3Di technology

Ankle fracture LP system

## Anatomic syndesmotic fixation

Fibula plates feature offset syndesmosis slots, allowing for centralized tibial fixation and anatomic reduction

## Comprehensive ankle fracture plating options


System features multiple plate designs that address variations in patient anatomies and fracture patterns

## Gravity™ Synchfix™

Syndesmosis fixation device

Medial button attachment technique preserves bone and reduces neurovascular impingement

Offset fibula plate

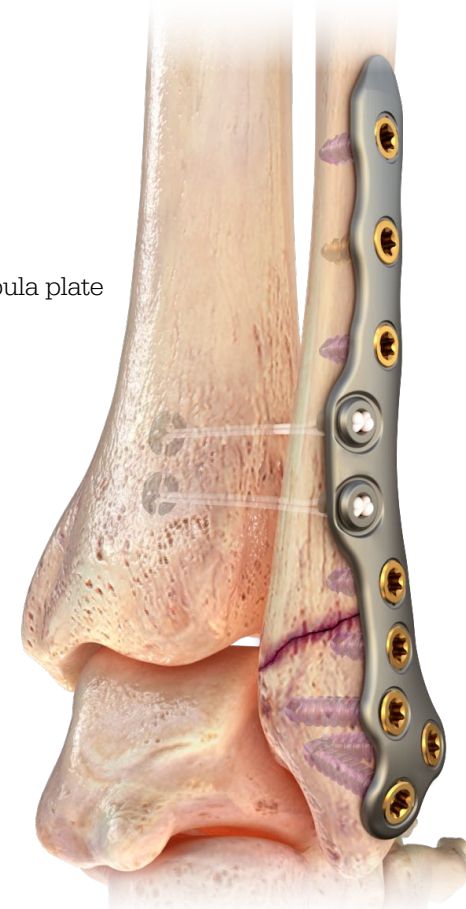

# Ortholoc 2

with 3Di technology

Crosscheck™ plating system

## Thinnest profile with an integrated lag screw<sup>1</sup>

The Ortholoc 2 Crosscheck is a diverse foot reconstruction system featuring the lowest profile plate with an integrated cross screw

- Predictable lag screw placement
- No bone removal (i.e. no counterbore)
- Fewer steps to improve procedural efficiency and save time in the O.R.

## Compression in lesser TMT fusions

The Y-plate is a utility solution with an integrated lag screw that delivers robust compression during lesser tarsometatarsal fusions



MTP plate



Y-plate

## AlloMatrix™ DBM putties

- Multiple formulations with and without cancellous bone chips
- Increased scaffolding for bone regeneration
- Moldable, packable bone graft substitute



# Ortholoc 2

with 3Di technology

Small bone plating system

## Diverse lesser metatarsal solutions

A broad range of utility and 5th metatarsal plate options provide versatility when faced with intraoperative challenges

## The market's first and only pinch plate<sup>2</sup>

Designed specifically for metatarsal neck fractures, the pinch plate's barbed arms grab and reduce smaller, hard-to-manage capital fragments



Hook plate


Pinch plate

## Ignite™ Power mix

- Powerful osteogenic graft (when mixed with BMA)
- Fracture callus potential
- Minimally invasive application



### Ankle fracture plating options

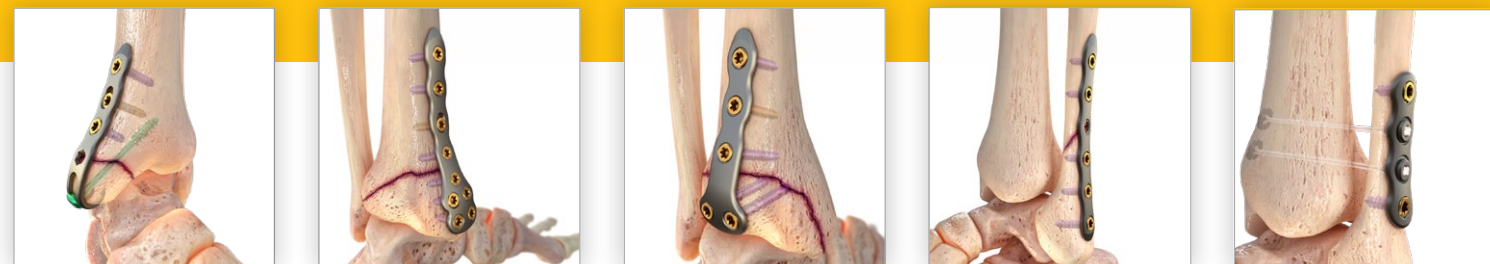


Hook plate    Medial tibial plate    Posterior tibial plate    1/3 tubular plate    Syndesmosis plate

### Crosscheck plating options



MTP plate    Lapidus 0mm    Lapidus 2mm    Utility plate 5 hole    Y-plate

### Small bone plating options

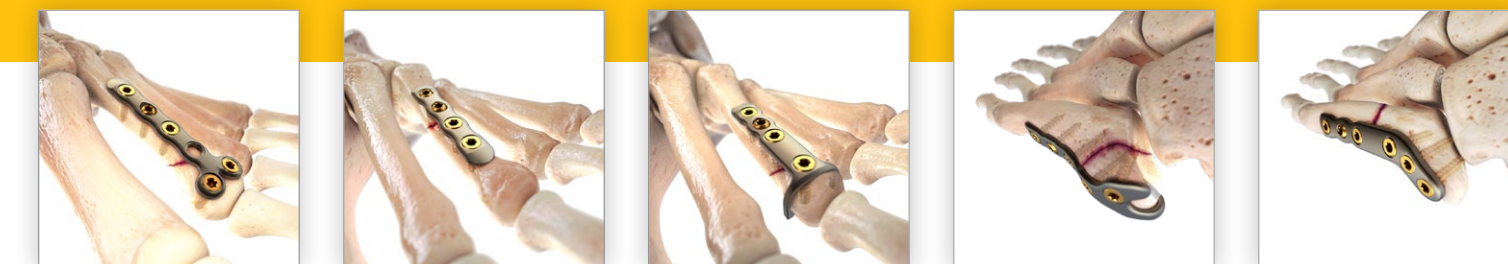

T-plate    Straight plate    Pinch plate    5th metatarsal Y-plate    Jones fracture plate

<sup>1</sup> Based on publicly available literature for Stryker, Paragon 28, Arthrex, Zimmer-Biomet, Depuy-Synthes, OsteoMed.

<sup>2</sup> In US market as of Feb 2018