


stryker


Driving
performance,

**delivering
versatility**

Spetzler-Malis and SILVERGlide Non-Stick Bipolar Forceps

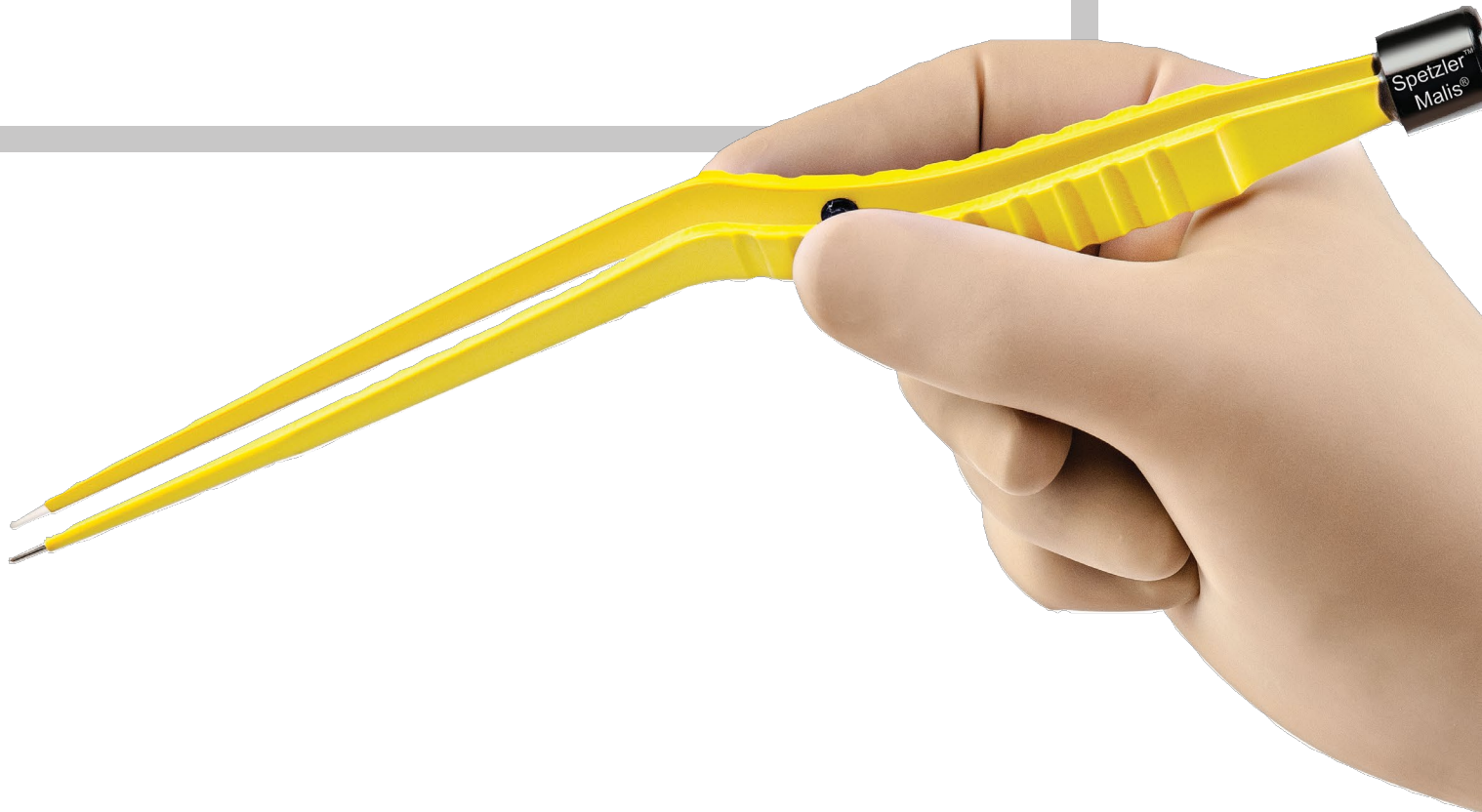

All of our non-stick bipolar forceps **resist tissue build-up** and are designed for **excellent visualisation**. Regardless of your preference for single-use Spetzler-Malis or reusable SILVERGlide, we can provide you with bipolar forceps engineered for **precision and control**.

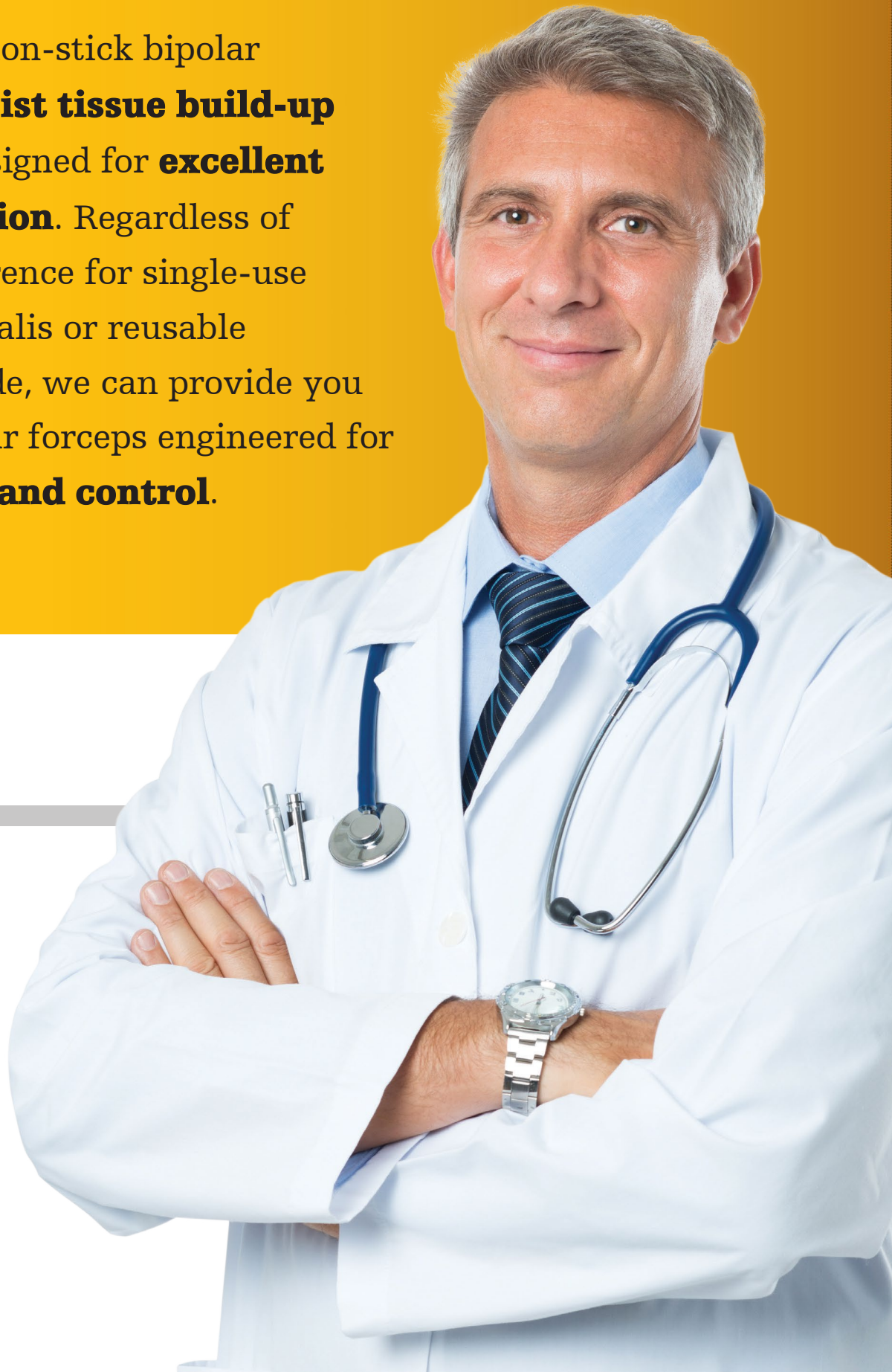


Achieve balance


and improve visualisation
with Spetzler-Malis Disposable
Non-Stick Bipolar Forceps

Disposable Spetzler-Malis Bipolar Forceps offer fresh advantages from every new use

- Single use ensures precise tip alignment and repeatable coagulation performance
- Non-stick technology enhances performance
- Sleek bayonet design improves visualisation
- Flexible lightweight cord and remote connection plug helps achieve balance
- Spring force elevates control and tactile feedback
- Transparent fine tubing of irrigating bipolar forceps allows for better observation
- Dual irrigation design of irrigating bipolar forceps maintains moisture of the tips and vessels

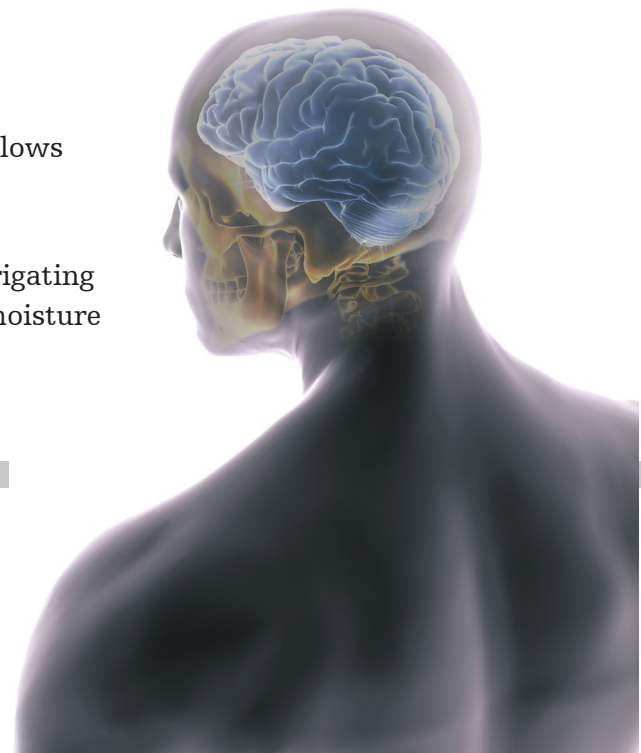

Disposable Spetzler-Malis **Standard** Bipolar Forceps

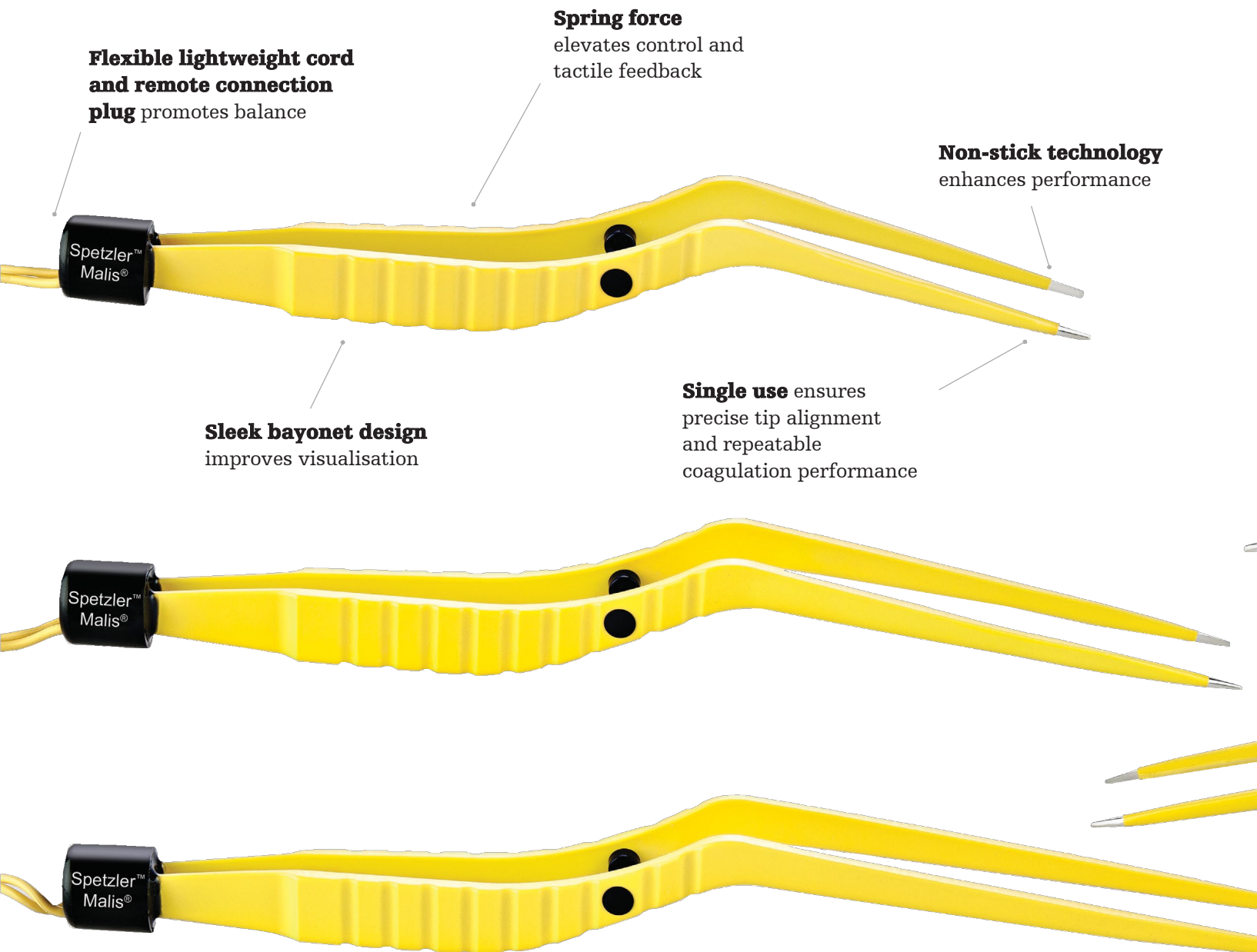

Flexible lightweight cord and remote connection plug promotes balance

Spring force elevates control and tactile feedback

Non-stick technology enhances performance

Sleek bayonet design improves visualisation

Single use ensures precise tip alignment and repeatable coagulation performance

18cm

6760-180-015 1.5mm tip
6760-180-010 1.0mm tip
6760-180-005 0.5mm tip

20cm

6760-200-015 1.5mm tip
6760-200-010 1.0mm tip
6760-200-005 0.5mm tip

23cm

6760-230-015 1.5mm tip
6760-230-010 1.0mm tip
6760-230-005 0.5mm tip

Disposable Spetzler-Malis **Dual Irrigating** Bipolar Forceps

Dual irrigation design
maintains moisture of the
tips and vessels

**Transparent fine
tubing** allows for
better observation



18cm

6770-180-015 1.5mm tip
6770-180-010 1.0mm tip
6770-180-005 0.5mm tip


20cm

6770-200-015 1.5mm tip
6770-200-010 1.0mm tip
6770-200-005 0.5mm tip

23cm

6770-230-015 1.5mm tip
6770-230-010 1.0mm tip
6770-230-005 0.5mm tip

Disposable Spetzler-Malis **Slim** Bipolar Forceps

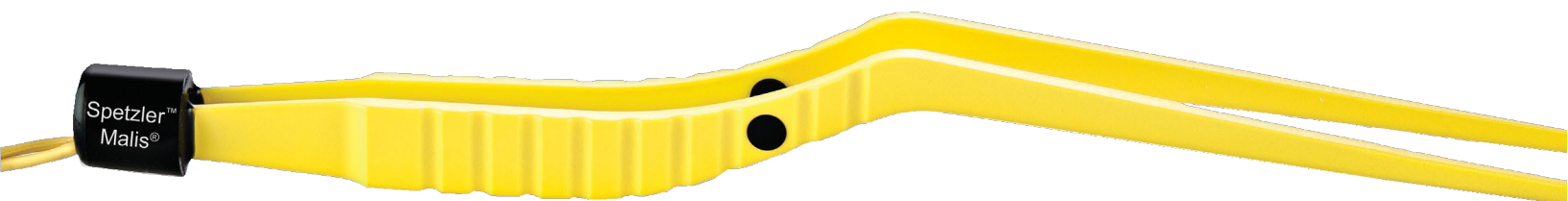


20cm

6780-200-015 1.5mm tip
6780-200-010 1.0mm tip
6780-200-005 0.5mm tip

23cm

6780-230-015 1.5mm tip
6780-230-010 1.0mm tip
6780-230-005 0.5mm tip

Slim profile provides access in narrow fields

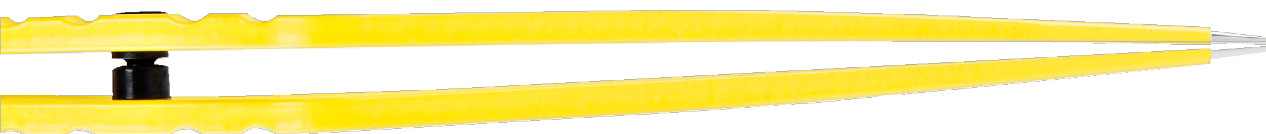

Disposable Spetzler-Malis Bipolar Forceps cord and irrigation sets

Rotary cord and tubing sets

6790-100-004 Bipolar irrigator integrated tubing and cord set, disposable, unitised plug

Gravity style cord and tubing sets

6790-100-002 Bipolar irrigator integrated tubing and bipolar cord set, disposable, unitised plug

Cord only sets

6790-100-006 Bipolar cord set, disposable, unitised plug




**Disposable Spetzler-Malis
Bipolar Forceps offer
fresh advantages from
every new use**

.5mm

1.0mm

1.5mm

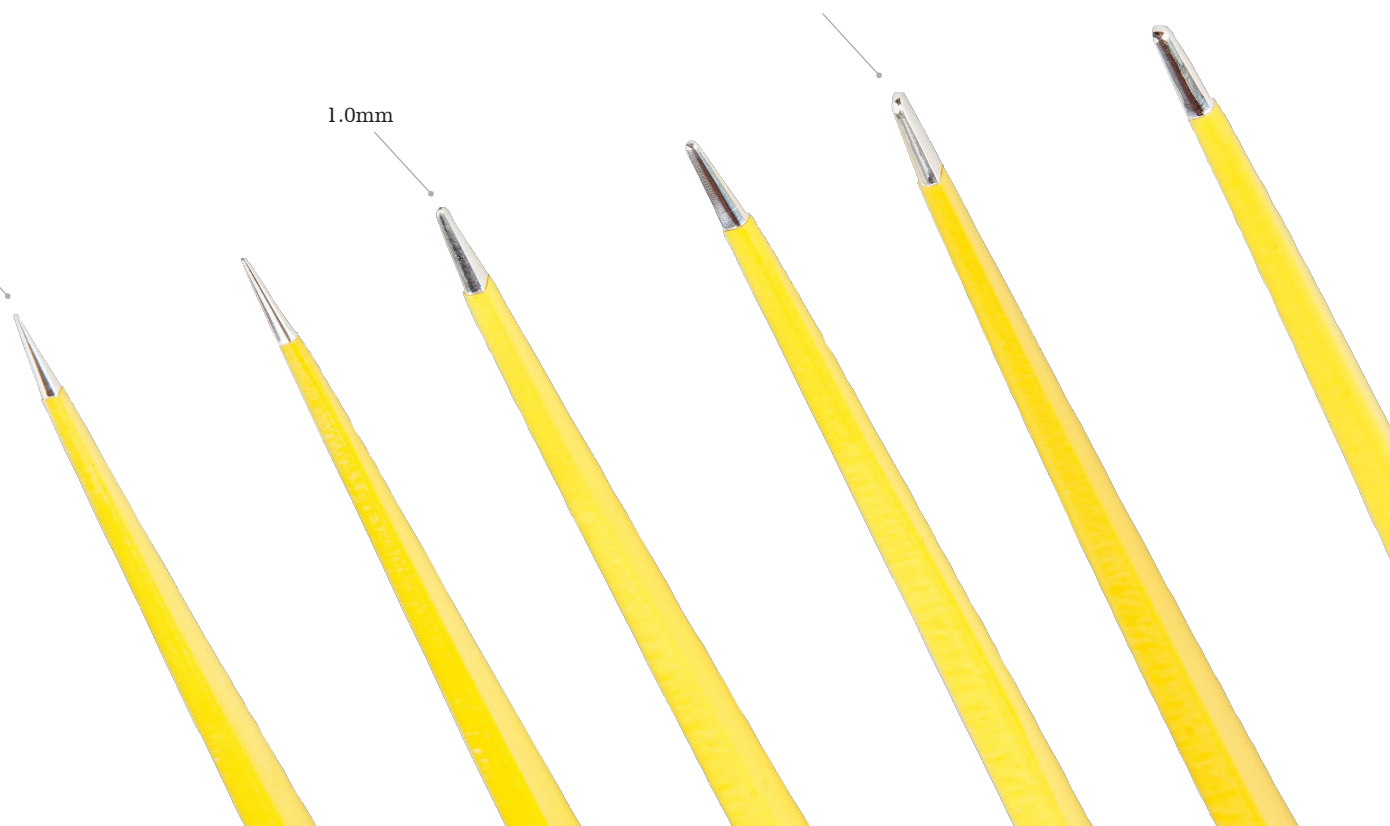

SILVERGlide Non-Stick ProSeries Bipolar Forceps



14.5cm bayonet

6730-145-015 1.5mm
6730-145-012 1.2mm
6730-145-010 1.0mm
6730-145-007 0.7mm
6730-145-004 0.4mm

16cm bayonet

6730-160-012 1.2mm
6730-160-010 1.0mm
6730-160-007 0.7mm
6730-160-004 0.4mm
6730-160-504, 45° up 0.4mm

19cm bayonet

6730-190-020 2.0mm
6730-190-015 1.5mm
6730-190-012 1.2mm
6730-190-010 1.0mm
6730-190-007 0.7mm
6730-190-004 0.4mm
6730-190-512, 45° up 1.2mm
6730-190-612, 45° down 1.2mm
6730-190-504, 45° up 0.4mm

22cm bayonet

6730-220-020 2.0mm
6730-220-015 1.5mm
6730-220-012 1.2mm
6730-220-007 0.7mm
6730-220-004 0.4mm
6730-220-512, 45° up 1.2mm
6730-220-612, 45° down 1.2mm
6730-220-507, 45° up 0.7mm
6730-220-607, 45° down 0.7mm
6730-220-307, 30° up 0.7mm
6730-220-107, 15° up 0.7mm

24cm bayonet

6730-240-020 2.0mm
6730-240-015 1.5mm
6730-240-012 1.2mm
6730-240-007 0.7mm
6730-240-004 0.4mm
6730-240-807, 80° up 0.7mm
6730-240-907, 80° down 0.7mm
6730-240-507, 45° up 0.7mm
6730-240-607, 45° down 0.7mm

28cm bayonet

6730-280-015 1.5mm
6730-280-012 1.2mm

19cm straight bayonet

6735-190-010 1.0mm
6735-190-007 0.7mm
6735-190-320, 30° up 2.0mm
6735-190-310, 30° up 1.0mm



Silver plating reduces collateral heating, "hot spots" and sticking for efficient performance

Ultra-thin tines allow better access in surgery





Elevating visibility through a slender design







Standard two-pin round connector works with a wide range of bipolar electrosurgical cables and generators

Lightweight handle for easy manipulation and balance

Guide-stop reduces scissoring and increases control

What sets SILVERglide apart?

The elimination of tissue build-up allows you to work without interruption. These non-stick characteristics continue to perform use after use. The SILVERglide alloy's level of thermal diffusivity is 37 times higher than that of stainless steel commonly used in forceps.

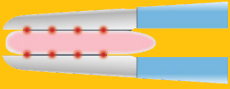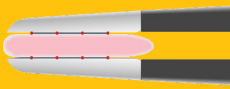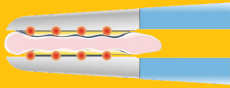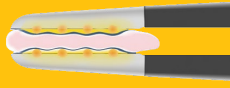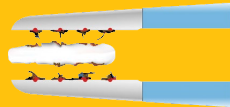
Standard stainless steel forceps	SILVERglide Non-Stick
<p>"Hot spots" from sparks occur during coagulation</p> 	<p>"Hot spots" start to occur during coagulation</p> 
<p>The "hot spots" cause localised sticking</p> 	<p>Heat is diffused rapidly to prevent sticking</p> 
<p>Tissue seared to the forceps tips is ripped away</p> 	<p>The forceps can be removed without tissue sticking</p> 

○ **Handle:** Light-weight, easy to grip and manipulate

○ **Guide-stop:** Reduces scissoring and increases control.

○ **Tines:** Ultra-thin tines allow better access in surgery.

○ **Tips:** Silver plating assists in dissipating heat quickly, reducing collateral heating and hot spots that cause sticking.

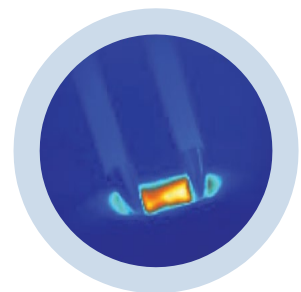

Malis Bipolar **Generator and Irrigator system**



Irrigator
Optional integrated irrigator

Malis Bipolar Generator and Irrigator system

- 6700-313-500 Malis Bipolar Generator
- 6700-313-600 Malis Bipolar Irrigator (does not include power cord)
- 6700-100-001 Bipolar Irrigator power cord
- 6700-313-004 Bipolar generator/irrigator interconnect cable
- 6700-313-002 Bipolar generator cut/coag dual foot pedal
- 6700-313-700 Bipolar generator cart

Other accessories

- 6700-313-003 Single foot pedal
- 6700-313-304 Single foot pedal adapter cable

Oversized colour screen
for easy comprehension

DUAL WAVE technology
provides confidence during cutting and coagulating

Programmable user profile stores information for up to eight surgeons

CUT and COAG knobs facilitate quick setting adjustments

Independent bipolar outputs enhance functionality



Whether you choose reusable or disposable, our non-stick bipolar forceps are designed to act as a specialised extension of your hands and eyes.

NSE

This document is intended solely for the use of healthcare professionals.

A healthcare professional must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that healthcare professionals be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A healthcare professional must always refer to the package insert, product label and/or instructions for use before using any Stryker product.

Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Malis, SILVERGlide, Spetzler-Malis and Stryker. All other trademarks are trademarks of their respective owners or holders.

The products depicted are CE marked in accordance with applicable EU Regulations and Directives.

For more information, please contact your local Stryker representative or visit Stryker.com

(9100004379EN)

9100004379EN Rev. A

DFS 09/2017

2017-13266

Copyright © 2017 Stryker

www.stryker.com



Stryker Instruments
4100 East Milham Avenue
Kalamazoo, MI 49001 USA
t: (+1) 269 323 7700
f: (+1) 800 999 3811
stryker.com



0197