

**Stryker I.D. Touch™ software:
A difference you can feel**

A surgeon's "touch" is one of his/her greatest assets, and our industry—exclusive I.D. Touch software aims to support this. Within the CORE™ Console power source, I.D. Touch software enables you to adjust torque from 0-100% to customize how your electric drill responds when pressure is applied during use—from low to high sensitivity—redefining performance to your liking. Your NSE sales representative can guide you through the options, and save the performance and feel that's right for you in user preferences.



Signature Portfolio Elite Attachments

Elite Attachments		
Part Number	Product Description	Color Band
5407-120-450	Elite 7cm Straight Attachment	Black
5407-120-452	Elite 7cm Angled Attachment	Black
5407-120-470	Elite 12cm Straight Attachment	Black
5407-120-472	Elite 12cm Angled Attachment	Black
5407-120-480	Elite 14cm Straight Attachment	Black
5407-120-482	Elite 14cm Angled Attachment	Black
5407-120-484	Elite 14cm Hooded Attachment	Black
5407-120-485	Elite 17cm Straight Attachment	Black
5407-120-486	Elite 17cm Angled Attachment	Black
5407-120-490	Elite 20cm Straight Attachment	Black
5407-120-492	Elite 20cm Angled Attachment	Black
Heavy Duty (HD) Attachments		
5407-120-050	HD 9cm Straight Attachment	Blue
5407-120-052	HD 9cm Angled Attachment	Blue
5407-120-070	HD 14cm Straight Attachment	Green
5407-120-072	HD 14cm Angled Attachment	Green
Irrigation Clips		
5407-010-270	Elite 7cm Irrigation Clip	N/A
5407-010-272	Elite 12-20cm Irrigation Clip	N/A
5407-010-200	HD 9cm Irrigation Clip	N/A
5100-010-200	Universal Irrigation Clip	N/A

To order or view the full line of Elite accessories and attachments, see your NSE sales representative, refer to your NSE product catalog, visit www.stryker.com/NSE or call **800-253-3210**.

NSE
Neuro Spine ENT

We're laser focused on Neuro, Spine and ENT. Our instruments, technology, education and support are all tailored to your specialties to help fulfill your unique needs and potentially make a positive difference in your day and outcomes.

- High Speed Drills (electric and pneumatic)
- Attachments and cutting accessories
- Sonopet® Ultrasonic Aspirator
- Bone Mill
- SilverGlide® Bipolar Forceps
- Microdebrider
- NasoPore® Bioresorbable Nasal Dressings
- HemoPore® Hemostatic Bioresorbable Nasal Dressing
- OtoPore® Bioresorbable Ear Dressings
- Surgeon headlights and other accessories in support of NSE specialties



Reconstructive

- Hips
- Knees
- Trauma & Extremities
- Joint Preservation
- Orthobiologics

Medical & Surgical

- Power Tools & Surgical Accessories
- Image Guided Navigation
- Endoscopy & Arthroscopy
- Integrated Communications
- Beds, Stretchers & EMS
- Sustainability Solutions

Neurotechnology & Spine

- Craniomaxillofacial
- Interventional Spine
- Neurosurgical, Spine & ENT
- Neurovascular
- Spinal Implants

The information presented in this brochure is intended to demonstrate the breadth of Stryker product offerings. Always refer to the package insert, product label and/or user instructions before using any Stryker product. Products may not be available in all markets. Product availability is subject to the regulatory or medical practices that govern individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Products referenced with TM designation are trademarks of Stryker. Products referenced with © designation are registered trademarks of Stryker.

Literature Number: 9100-003-173 Rev None G55/PS

Copyright © 2015 Stryker
Printed in USA

Stryker Neuro Spine ENT
4100 East Milham Avenue
Kalamazoo, MI 49001 USA
t: 269 323 7700 f: 800 999 3811
toll free: 800 253 3210

www.stryker.com/NSE



Neuro Spine ENT



Elite Attachments



The Elite line

Expanding our signature style for the Signature Portfolio

With its industry-unique telescoping design, you've made Elite attachments and accessories among our most popular. Now, your options have expanded even further with a refined line of attachments that deliver a more natural balance, grip and maneuverability.


As always, Elite instruments aim to provide the precise functionality you need while also streamlining ease-of-use and inventory. Elite cutting accessories are compatible with any Elite attachment—and can be adjusted up to five exposure levels (up to 10mm). This reduces the need to change attachments or cutting accessories to achieve the reach you want, and can also reduce inventory items up to 5:1.



5407-120-XXX

Patented multi-notch telescoping cutting tools enable convenient customization of exposure setting for specific user preferences and applications

- Adjustable up to five exposure levels (up to 10mm)
- Compatible with any Elite attachment
- Reduces and simplifies inventory
- Streamlines set up and packaging waste
- Reduces opening multiple cutting accessories for different attachments

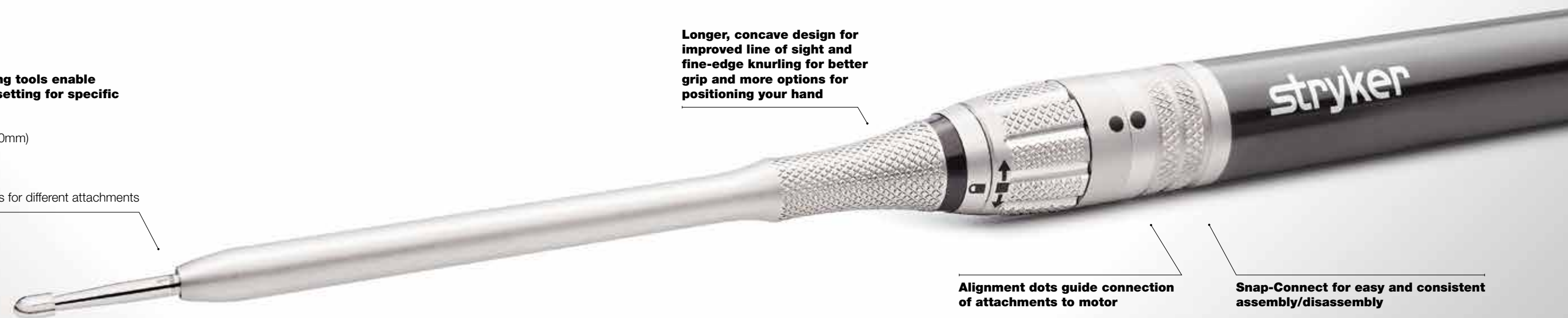

Longer, concave design for improved line of sight and fine-edge knurling for better grip and more options for positioning your hand

Alignment dots guide connection of attachments to motor

Snap-Connect for easy and consistent assembly/disassembly

Easily identifiable marks help surgeon and staff adjust exposure to desired level

Refined attachment design for improved balance and maneuverability

Portfolio options

Below are some of the style options of the nearly 200 Elite cutting accessories that can help you achieve the specific functionality you want. For the full line, see your NSE sales representative, NSE product catalog or materials on www.stryker.com/NSE.

Precision Match Head



Match Head Less Aggressive



Match Head Soft Touch



Precision Round



Match Head



Round Fluted Aggressive



Round Fluted Soft Touch



Round Diamond

