

Sage PrimaFit[®]

External Urine Management for the Female Anatomy

Easy connection

Suction tubing adapter allows for easy connection to hospital tubing

Accommodates a range of body types


Long, flexible tubing allows multiple suction tubing angles

Gentle adhesive pad

Keeps device in place and can be easily and gently removed

Helps keep skin dry

Ultra-soft fabric diverts urine away from skin

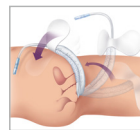


Flexible fit

Contours to patient anatomy and maintains shape for the duration of use

Stays in place

Soft, tapered end gently fits in perineal area

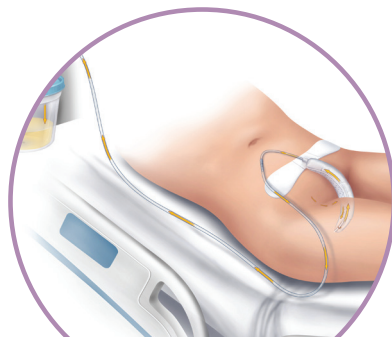

Helps promote
early catheter removal

Address the #1 CAUTI risk factor, the indwelling catheter¹

PrimaFit is an external urine management system designed to help manage urinary incontinence in patients with female anatomy and promote early catheter removal, addressing the #1 risk factor of catheter-associated urinary tract infections (CAUTI). PrimaFit serves as an alternative to incontinence pads, diapers, or, at times, catheters while transitioning patients to independent continence.

When properly positioned, PrimaFit's ultra-soft fabric diverts urine away from the patient's skin. Urine is then absorbed into the system's core and suctioned into a collection canister, helping keep the skin dry and addressing a risk factor or incontinence-associated dermatitis (IAD).

No catheter. No CAUTI.[™]

Approximately

75%

of hospital-acquired UTIs are associated with an indwelling urinary catheter.¹



(1) PrimaFit, adult size/package

40 devices/case
Reorder #5400