

Cayman[®]

United Plate System




The Cayman United Plate System is specifically designed for rigid fixation to the Cascadia Lateral 3D line of interbody implants. The single level, one- and two-hole plate configurations allow for a customizable construct from the lateral approach to help achieve sagittal balance. The Cayman United plate utilizes the revolutionary SMARTLock Technology allowing bone screws to lock directly to the plate without an additional locking mechanism.

Cayman United Plate System


Single-level plate

- **Single-level plate sizes** are designed to minimize retraction required for placement
- **Assembly alignment feature** keys into Cascadia interbody providing resistance to rotation
- **One- and two-hole plate configurations** sized to all Cascadia Lateral height and lordosis options

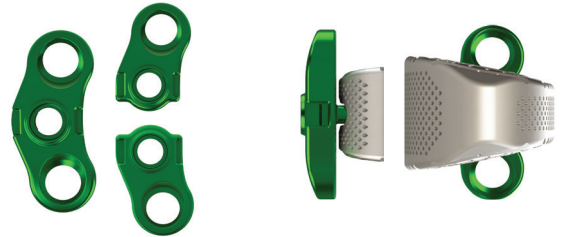

Versatile instrumentation

- **Screw-thru Plate Inserter** allows for either all-in-one insertion or in-situ assembly after the interbody has been placed
- **15° nominal screw angulation** with SMARTLock technology allows for screw placement parallel with endplates to 30° offset, allowing freedom to work within a limited access corridor

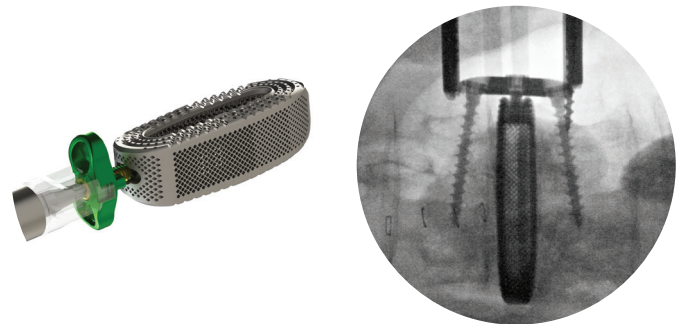

Additional features

- **Complete self-drilling screw offering** available in both 5.0 and 5.5mm diameters, with lengths of 28, and 40 – 60mm (in 4mm increments)
- **One-step locking** using the 20 in-lb Torque Limiting Handle for both bone screws and plate assembly screws



Spine division

A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A surgeon must always refer to the package insert, product label and/or instructions for use before using any Stryker product. The products depicted are CE marked according to the Medical Device Directive 93/42/EEC. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Cascadia, Cayman, K2M, SMARTLock, Stryker. All other trademarks are trademarks of their respective owners or holders.

CAYUN-SS-2_23126
Copyright © 2020 Stryker



K2M, Inc.
600 Hope Parkway SE
Leesburg, VA 20175
t: 571 919 2000

www.stryker.com