

Chesapeake[®]


Cervical-Ti Stabilization System



The Chesapeake Cervical-Ti Stabilization System is a multi-screw construct designed to provide stability to the anterior column while reducing the need for supplemental fixation in the cervical spine. Featuring SMARTLock Technology, the system exemplifies advances in new and innovative instrumentation.


Chesapeake Cervical-Ti Stabilization System

Implant sizes

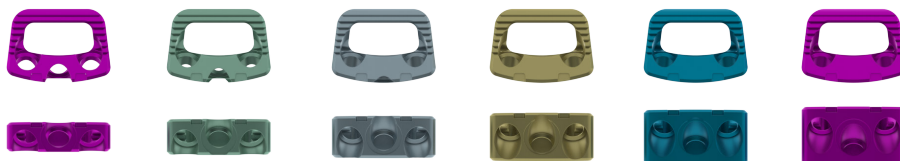

5 – 10mm heights in 1mm increments. Color-coded for ease of identification.


Screw lengths - Ø3.5mm and Ø3.85mm

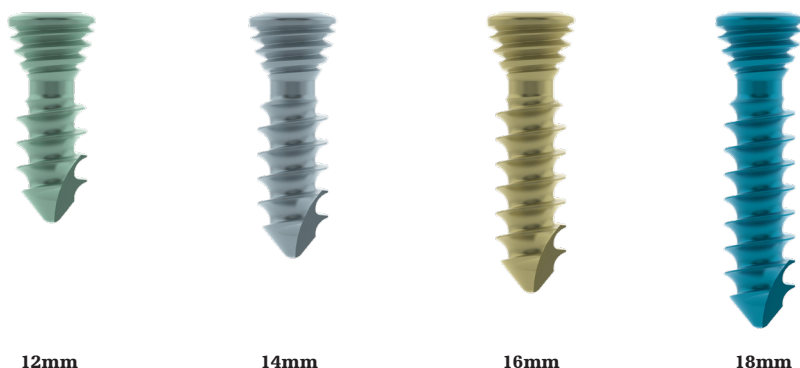

Instrument design features



- All-in-one Drill Guide Inserter is used to awl, drill, tap, and insert screws at a fixed trajectory of 40°
- Adjustable stop inserter provides the ability to adjust the depth of the interbody into the intervertebral disc space
- Awls, Drills, Taps, and Screw drivers are available in both fixed and angled options
- Color-coded to match the implant size for ease of identification intraoperatively

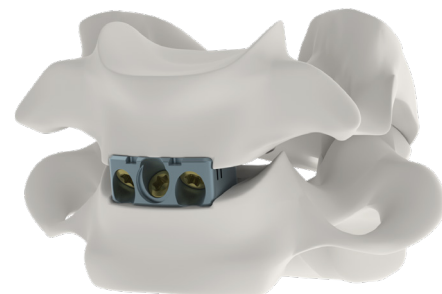

Chesapeake Cervical-Ti
Stabilization System

- Featuring SMARTLock Technology, which requires no additional locking mechanism
- 12 x 17 and 13 x 18mm footprints
Available in 7° lordotic options
- Anterior fixation with zero-profile design
- A multi-screw construct with 14° converging angles
- Interbodies manufactured from CP Grade II titanium
- Screws manufactured from titanium alloy (Ti6Al4V)

Spine division

A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A surgeon must always refer to the package insert, product label and/or instructions for use before using any Stryker product. The products depicted are CE marked according to the Medical Device Directive 93/42/EEC. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Chesapeake, K2M, SMARTLock, Stryker. All other trademarks are trademarks of their respective owners or holders.

CHECTI-SS-1_23809
Copyright © 2020 Stryker



K2M, Inc.
600 Hope Parkway SE
Leesburg, VA 20175
t: 571 919 2000

www.stryker.com