


# Ozark<sup>®</sup>

## Guide Cervical Plate System



The Ozark Guide Cervical Plate System features a per-level, integrated locking cover that provides visual and tactile confirmation of the final lock position. The Guide Plate features a fixed All-In-One Drill Guide through which the surgeon can awl, drill, tap and place screws. The system also includes a variety of plate and screw sizes, and surgical instrumentation to create constrained, semi-constrained, or hybrid screw constructs.

# Ozark Guide Cervical Plate System

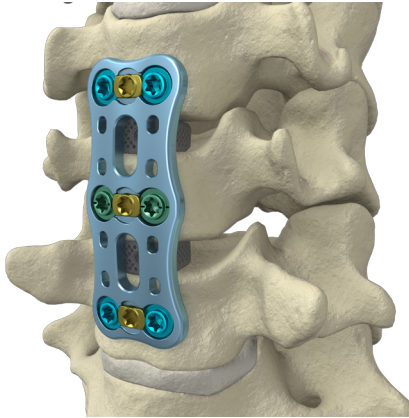


## Plate design

- 1–5 level plates with integrated locking covers
- Clear visual and tactile confirmation of final lock position
- 17.4mm wide and 2.5mm thin
- Pre-lordosed to help match patient anatomy
- Strategically placed windows for graft visualization



Fixed screws

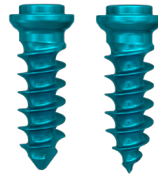

Variable screws

## Screw design

- Self-starting and self-tapping screw options
- Color-coded by fixed and variable, as well as by screw length
- Ø4.0mm standard and Ø4.5mm rescue screws
- 12, 14, and 16mm standard lengths (additional lengths available by request)
- Flexibility to insert variable screws up to 30° conically



Variable drill guide



Fixed drill guide



FAST guide

## Instruments

- Same driver used for screw insertion and final locking
- Fixed angle screw-thru (FAST) guide provides drill, awl, tap, and screw-thru option
- Multiple drill guide lengths to accommodate surgeon preference

## Spine division

A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A surgeon must always refer to the package insert, product label and/or instructions for use before using any Stryker product. The products depicted are CE marked according to the Medical Device Directive 93/42/EEC. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: K2M, Ozark, Stryker. All other trademarks are trademarks of their respective owners or holders.

OZRGD-SS-1\_25305  
Copyright © 2021 Stryker



K2M, Inc.  
600 Hope Parkway SE  
Leesburg, VA 20175  
t: 571 919 2000

[www.stryker.com](http://www.stryker.com)