

SharpShooter

**Inside-out meniscal
repair system**



**23 gauge needle and
2-0 high strength suture**



**360-degree cannula rotation
allows for customizable
hand positions**

**Simple snap mechanism for
loading cannulas**



**Anatomic
cannulas**

**Trigger delivery
is designed to
incrementally advance
needles single-handedly**


**Sterile,
single-use
disposable
handle**

SharpShooter

**Stryker's inside-out
meniscal repair system**

A simple, easy-to-use method for repairing the meniscus utilizing an inside-out suture technique. The system includes a unique method to deliver either high strength 2-0 UHMWPE suture or 2-0 polyester braided suture using a delivery handle and six anatomic cannulas to reach all areas of the meniscus arthroscopically. The SharpShooter is designed to give enhanced control of needle placement and delivery.

Popliteal retractor
**This specially
designed shield
deflects needles away
from neuro-vascular
structures during
inside-out repairs**

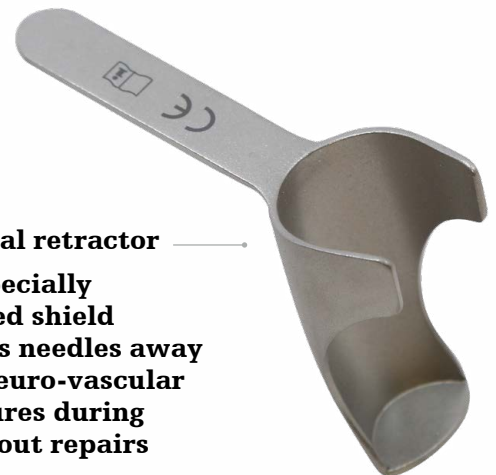

Features and benefits

Targeted delivery

- SharpShooter handle—A disposable suture passer that incrementally advances the SharpShooter needle
- SharpShooter cannulas—Six specially designed cannulas easily attach to the SharpShooter handle allowing access to all zones of the meniscus
- Designed to provide complete control of targeting needles

High-strength suture

- White 2-0 UHMWPE and green 2-0 polyester braided suture provide color options while maintaining the necessary strength

Efficient design

- Single-handed needle advancement requires less operating assistance
- Double-armed, trocar-tipped suture needle incorporates visual markings designed to distinguish the distance traveled through the cannula before engaging the meniscus



SharpShooter offers angles to reach many locations

Part Number	Description
4700	SharpShooter handle
PB4704	SharpShooter surgical suture, braided PE
HS4704	SharpShooter surgical suture, UHMWPE
4750	SharpShooter cannula right posterior
4751	SharpShooter cannula left posterior
4752	SharpShooter cannula right middle
4753	SharpShooter cannula left middle
4754	SharpShooter cannula right anterior
4755	SharpShooter cannula left anterior
233050115	Popliteal retractor small
233050116	Popliteal retractor large

Sports Medicine

This document is intended solely for the use of healthcare professionals. A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. We do not dispense medical advice and recommend that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate Stryker's products. A surgeon must always refer to the package insert, product label and/or instructions for use, including the instructions for cleaning and sterilization (if applicable), before using any of Stryker's products. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets.

Please contact your representative if you have questions about the availability of Stryker's products in your area.

Stryker or its affiliated entities own, use, or have applied for the following trademarks or service marks: SharpShooter and Stryker. All other trademarks are trademarks of their respective owners or holders.

The absence of a product, feature, or service name, or logo from this list does not constitute a waiver of Stryker's trademark or other intellectual property rights concerning that name or logo.

1000902220 Rev C

Copyright © 2019 Stryker

Stryker
Sports Medicine
5670 Greenwood Plaza Blvd.
Ste. 200
Greenwood Village, CO 80111
t: 866 596 2022
www.sportsmedicine.stryker.com