

Evolve® EPS

Elbow Plating System (EPS)



Medial/Posterolateral – 90°
(orthogonal)



Medial/Lateral – 180°
(parallel)



Olecranon

Evolve EPS

Elbow Plating System

- Polished stainless steel
- Ortholoc polyaxial locking
- Anatomically contoured low profile plate design
- 90 or 180 degree humeral plating options

For more information, visit stryker.com



Evolve EPS

Elbow Plating System (EPS)

Anatomically contoured plates

- Low profile design
- Polished stainless steel material
- Option to bend if necessary

Comprehensive plate and screw options


- Orthogonal (90°) or parallel (180°) humeral plating options
- Multiple plate lengths
- 2.7mm cortical and 3.5mm cancellous and cortical screws compatible with all holes

Ortholoc polyaxial locking


- Allows for screw placement in 30 degree cone - 15 degrees from center axis
- No need to tap threads in plate intraoperatively

Enhanced instrumentation


- Instrumentation designed for ease of use
- **In situ** bending to more appropriately fit the bone
- Comprehensive instrument system in 2 trays

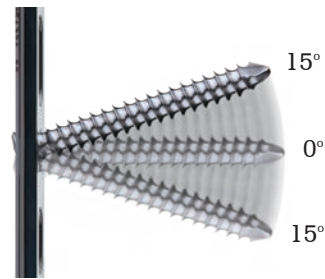

In Situ Bending

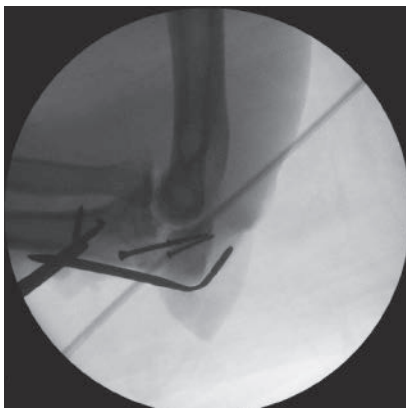

Pre bend

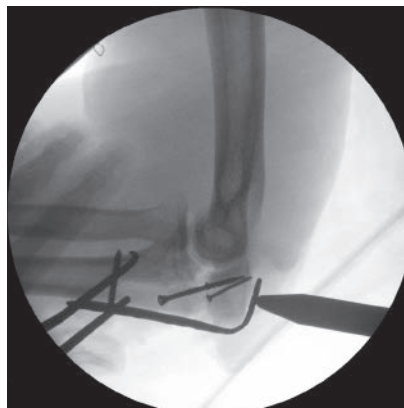

Bending



Post bend

Images courtesy
Dr. John Capo

This document is intended solely for the use of healthcare professionals. A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate a Stryker product. A surgeon must always refer to the package insert, product label and/or instructions for use, including the instructions for cleaning and sterilization (if applicable), before using any Stryker product. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Evolve, Stryker.

All other trademarks are trademarks of their respective owners or holders.

010440B_05-03-2021

Copyright © 2021 Stryker

Manufactured by:

Wright Medical Technology, Inc.
1023 Cherry Road
Memphis, TN 38117

161 Rue Lavoisier
38330 Montbonnot Saint Martin
France

stryker.com